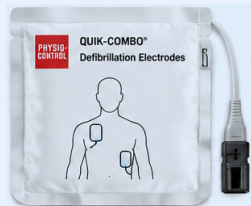
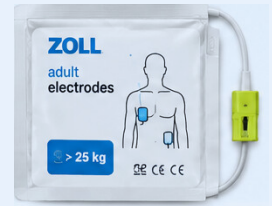
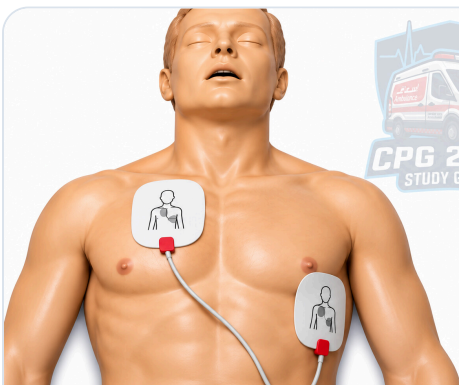


**1 Age → Mode**
Paediatric ≤13 y · Adult >13 y**2 Weight → Pads**
Device-specific kg cut-off**3 Device → Joules**
Never mix device doses**1 Mode Selection** age-based (UAE classification)**Paediatric****≤ 13 years****MANUAL MODE ONLY**

⚠ Do not press ANALYZE. AED mode follows the adult protocol (LIFEPAK 15: 200–300–360 J) — too much for children. If ANALYZE was pressed and the rhythm is shockable, select energy manually.

Adult**> 13 years****MANUAL OR AED MODE**

Fixed device doses — **LIFEPAK 15:** 200 > 300 > 360 J
ZOLL X-Series: 120 > 150 > 200 J

2 Pad Selection weight-based**LIFEPAK 15****Paediatric pads**
≤ 15 kg**Adult pads**
> 15 kg**ZOLL X-Series****Paediatric pads**
≤ 25 kg
or ≤ 8 years**Adult pads**
> 25 kg**3 Pad Placement** ensure pads adhere and never touch each other**Adult — anterolateral**

Right infraclavicular · left mid-axillary (apex)

**FRONT****BACK****Infant / small child — antero-posterior**

Centre of chest + between the shoulder blades, to avoid pad overlap

⚠ No paediatric pads? Do not delay defibrillation. Use adult pads and select the paediatric (weight-based) joules manually — justify in the ePCR why paediatric pads were not used. Never use paediatric pads above the device's weight limit.



4 Energy VF / pulseless VT · joules

LIFEPAK 15

ADULT · FIXED DOSE

1ST **200 J** > 2ND **300 J** > 3RD + **360 J**

PAEDIATRIC · WEIGHT-BASED, NEVER ABOVE

1ST **2 J/kg**
max 200

2ND **4 J/kg**
max 300

3RD + **4–10 J/kg**
max 360

ZOLL X-Series

ADULT · FIXED DOSE

1ST **120 J** > 2ND **150 J** > 3RD + **200 J**

PAEDIATRIC · WEIGHT-BASED, NEVER ABOVE

1ST **2 J/kg**
max 120

2ND **4 J/kg**
max 150

3RD + **4–10 J/kg**
max 200

5 Paediatric Ready Reckoner manual mode · 2 → 4 → 4–10 J/kg, capped per device

LIFEPAK 15

KG	PADS	1ST J	2ND J	3RD + J
3	Paed	6	12	12–30
4	Paed	8	16	16–40
5	Paed	10	20	20–50
6	Paed	12	24	24–60
7	Paed	14	28	28–70
8	Paed	16	32	32–80
9	Paed	18	36	36–90
10	Paed	20	40	40–100
12	Paed	24	48	48–120
14	Paed	28	56	56–140
15	Paed	30	60	60–150
16	Adult	32	64	64–160
18	Adult	36	72	72–180
20	Adult	40	80	80–200
22	Adult	44	88	88–220
25	Adult	50	100	100–250
28	Adult	56	112	112–280
30	Adult	60	120	120–300
35	Adult	70	140	140–350
40	Adult	80	160	160–360*
50	Adult	100	200	200–360*

ZOLL X-Series

KG	PADS	1ST J	2ND J	3RD + J
3	Paed	6	12	12–30
4	Paed	8	16	16–40
5	Paed	10	20	20–50
6	Paed	12	24	24–60
7	Paed	14	28	28–70
8	Paed	16	32	32–80
9	Paed	18	36	36–90
10	Paed	20	40	40–100
12	Paed	24	48	48–120
14	Paed	28	56	56–140
15	Paed	30	60	60–150
16	Paed	32	64	64–160
18	Paed	36	72	72–180
20	Paed	40	80	80–200
22	Paed	44	88	88–200*
25	Paed	50	100	100–200*
28	Adult	56	112	112–200*
30	Adult	60	120	120–200*
35	Adult	70	140	140–200*
40	Adult	80	150*	160–200*
50	Adult	100	150*	200*

200* capped at the device's adult dose for that shock 3rd + : start ≥4 J/kg, escalate up to 10 J/kg Bold line = pad change (15 kg LIFEPAK · 25 kg ZOLL)

Values rounded to the nearest joule

6 Operational Reminders

No ECG through the pads?

Press the LEAD button to select pads on both devices.

Manual shock sequence

Select energy > CHARGE > state "all clear" and confirm > SHOCK.

Document in ePCR

Age, pad type, mode (manual/AED), joules delivered and rhythm.

Monitored KPIs — paediatric arrest

Time to first shock · correct joules and adherence to the paediatric sequence.